

## Pacific Examination and Treatment Tables

FIJI P5, P7  
HAWAII S2, S4



This is medical device. Use it according to the  
user manual or label.  
The information contained in this leaflet is intended  
for healthcare professionals only.

**Pacific Examination and Treatment Table** is intended for standard examinations and procedures carried out in a wide range, also used for simple medical activities. It is usually used in professional medical care facilities, i.e. hospitals, clinics, in sanatoriums, nursing homes, as well as in doctor's offices, treatment rooms. Well-thought-out design, the possibility to adjust the angles of the sections and the stability of the product ensure that the patient is in a comfortable position for examination/treatment.

The benefits of treatment tables are based on the fact that it is possible for the patient to take the right position for the treatment, the possibility to adjust the height and/or sections of the table, which facilitates the work of the rehabilitator/physiotherapist and increases the comfort of the patient during the massage treatment. PACIFIC table is characterized by stability, reliability and ease of use. These features result from careful design, manufacturing methods applied in production and the quality assurance system.

### Pacific HAWAII features:

- 2- or 4-section table with electric height adjustment,
- hand-held remote control,
- head section adjustment by the means of gas spring, with shaped face hole with crescent plug,
- foldable armrests lowered till 19 cm (only in 4-section table),
- multi-layered and very soft upholstery of 15 cm thickness,
- Trendelenburg/reverse Trendelenburg position (15°),
- 2 wheels, 2 feet.



Pacific HAWAII

### Pacific FIJI features:

- 5- or 7-section table with electric height adjustment,
- hand-held remote control,
- head section adjustment by the means of gas spring, with shaped face hole with crescent plug,
- foldable armrests lowered till 19 cm (only in 7-section table),
- adjustment of: back / thighs / shank sections by means of electric actuators,
- multi-layered and very soft upholstery of 15 cm thickness,
- Trendelenburg/reverse Trendelenburg position (15°),
- 2 wheels, 2 feet.



Pacific FIJI

### Indications:

- upper crossed syndrome,
- delayed onset muscle soreness,
- neck pain,
- occupational overloads,
- muscular tension disorders,
- spinal disc disorder,
- osteoarthritis of the hip, knee joints,
- chronic lumbosacral pain,
- shoulder pain syndromes,
- cystic fibrosis,
- chronic back pain,
- chronic low back pain,
- pain in the quadriceps fatigue,
- tension-type headache.





#### Additional options:

- integrated option of heating and LED lights,
- central lift (2 castors central lifted).

# Tables articulations:



PACIFIC Table

type **FIJI**, model **P5**

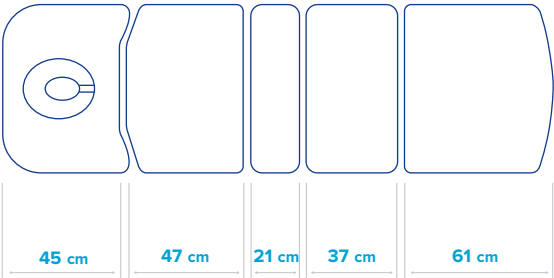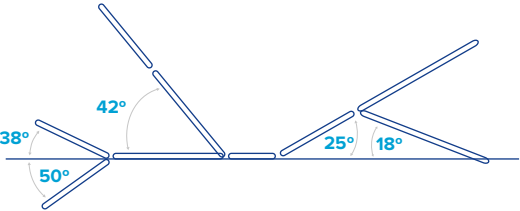

PACIFIC Table

type **FIJI**, model **P7**

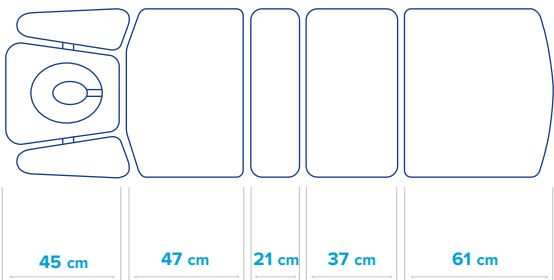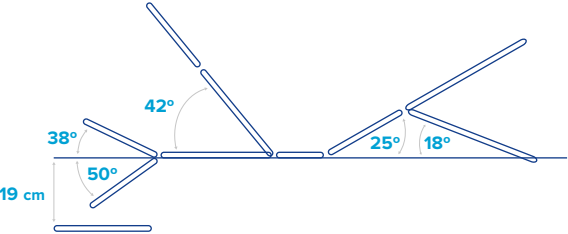


PACIFIC Table

type **HAWAII**, model **S2**




PACIFIC Table

type **HAWAII**, model **S4**

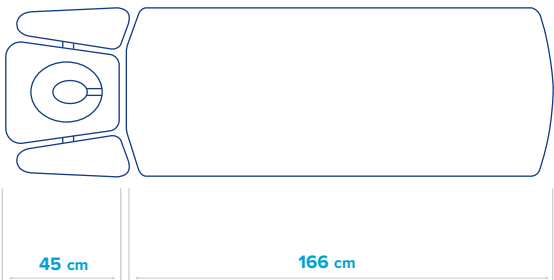

## Technical data:

Table parameters		PACIFIC			
		FIJI		HAWAII	
		P5	P7	S2	S4
Table length [cm]	total	212			
	head section	45			
	back section	47	166		
	middle section	21			
	thighs section	37			
	shank section	61			
Angle adjustment [°]: raising/lowering	head section	38:50			
	back section	42:00	-		
	thighs section	25:00	-		
	shank section	25:18	-		
Upholstery thickness	[cm]	15			
Angle adjustment: Trendelenburg and reverse Trendelenburg	[°]	+15:-15			
Height adjustment	[cm]	min. 64 ±1 max. 104 ±1			
Table width	[cm]	approx. 80			
Lifting capacity [kg]		≤ 150			
		(maximum safe load while adjusting height)			
Weight [kg]	not more than	150	130		
Power and security requirements	voltage, frequency	120-240V~ / 50/60 Hz			
	current consumption	max. 4,4 A			
	protection class	II			
	applied part	type B			



Central lift - option



Heating and LED lights - option



FIJI remote control