


Bathaid & Hydraid



Care bathtub **Bathaid** / Bathtub hoist **Hydraid**

This is medical device. Use it according to the
user manual or label.

The information contained in this leaflet is intended for
healthcare professionals only.

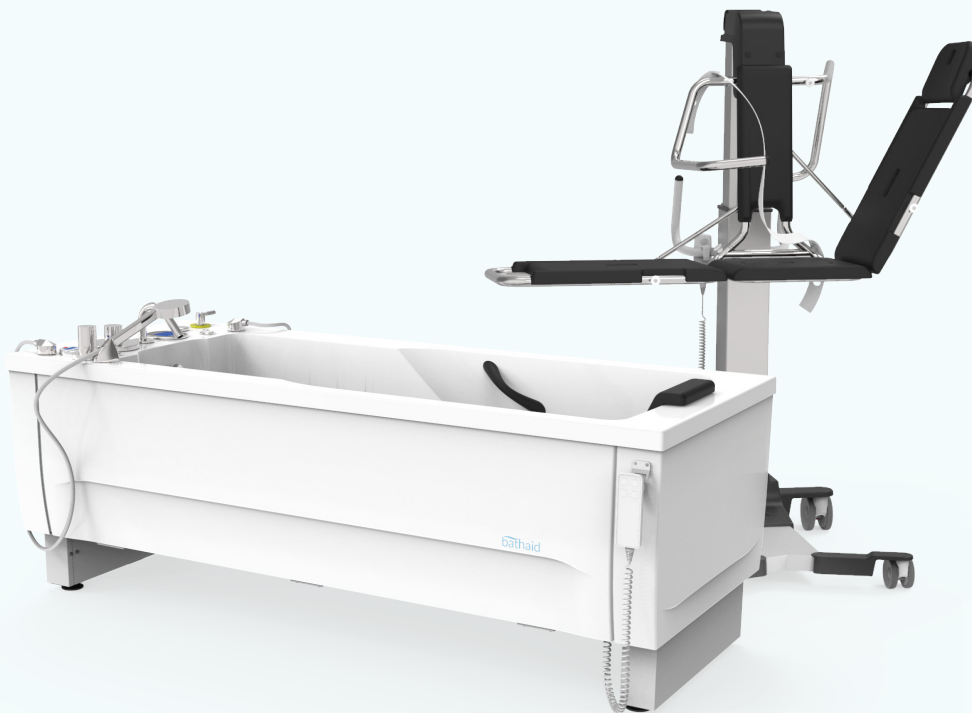

The care bathtub **Bathaid** is a modern, advanced design tub, specially created and intended for bathing and washing patients in hospitals and nursing homes, as well as in long-term care facilities. The built-in height adjustment mechanism makes it easier for the personnel to perform hygienic activities related to the daily care of the patient.

The Bathaid makes it easier for personnel to perform a nursing bath for a patient with reduced mobility due to injury, illness or disability.

Bathaid (basic equipment):

- ergonomic design,
- specially profiled basin with arm supports,
- automatic siphon,
- electric lifting and lowering of the bathtub basin with adjustment in the range of 40 cm,
- lifting and lowering the basin by the means of a fixed panel and a remote control (with access authorization), which allows the basin to be lifted and lowered with two different lifting/lowering speeds,
- built - in bathtub disinfection system,
- shower head for washing the patient,
- handles,
- headrest,
- leg support.

The bathtub hoist **Hydraid**, with height adjustment, is designed for use in hospitals, sanatoriums and nursing homes. It is used for lifting, transferring and bathing person with reduced mobility due to illness or dysfunction. It allows the caregiver to transfer them from the bed directly to the bathtub and back again.

It can be operated by one person, which allows for ergonomic work of the personnel.

The Hydraid can be used together with the care bathtub Bathaid.

Hydraid:

- raising/lowering,
- mobile,
- equipped with a control panel and a manual remote control for height adjustment,
- 3 - piece lying (+ pillow),
- adjustment of the lying sections.

Care bathtub **Bathaid**

2 models of the Bathaid:

ECONOMY



Additional equipment:

- bathtub spout with thermostat,
- electronic basin water temperature indicator.

Optional equipment:

- AUDIO

PREMIUM



Additional equipment:

- touch control panel with colour display,
- thermostatic tub faucet with an independent spout,
- automatic filling of the bathtub basin (one treatment level),
- electronic indicator of the water temperature in basin and spout.

Optional equipment:


- AUDIO,
- CHROMO - 12 LED points with selecting of lighting effects.

Bathtub hoist **Hydraid**

Hydraid (as seat)



+ option

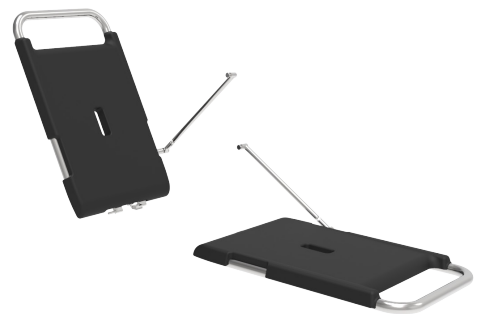


Hydraid (as lying)



Care bathtub **Bathaid**

Technical data:

Dimensions (L x W) [mm]:	1990 x 800
Height min./max (basin edge) [mm]:	720-1120
Basin depth [mm]:	430
Drain Ø [mm]:	50
Capacity up to overflow [l]:	270
Power supply [V, Hz]:	230, 50
Device weight (ECONOMY/PREMIUM) [kg]:	210/230
Max. patient weight [kg]:	185
Max. safe working load of the bathtub (SWL) [kg]:	400
Max. filling time up to overflow (ECONOMY/PREMIUM) (at flow ~70l/min) [min]:	~16/~9,5
Max. drainage time of the bathtub basin [min]:	~4
Lifting/lowering time (at SWL 400 kg) [s]:	29/23

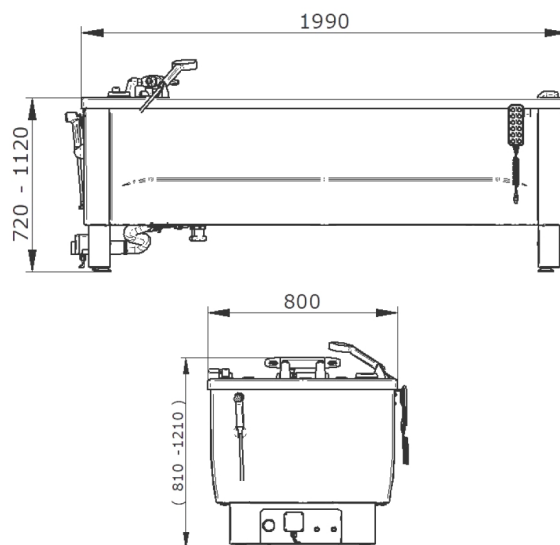

Bathtub hoist **Hydraid**

Technical data:

Dimensions (as seat) (L x W x H (min./max.)) [mm]:	970 x 760 x 1240/1840
Dimensions (as lying) (L x W x H (min./max.)) [mm]:	970 x 1840 x 1240/1840
Dimensions of lying top (L x W) [mm]:	430 x 1840
Seat width [mm]:	435
Height of grip [mm]:	1190
Turn diameter [mm]:	1150
Wheels diameter (with/without brake)[mm]:	125/100
Input voltage [V, Hz]:	100-240, 50-60
Battery capacity [Ah]:	2,9
Device weight [kg]:	95
Max. safe working load (SWL) [kg]:	160
Lifting speed (raising/lowering) [s]:	42/35

