

# AQUAI

Dry hydromassage membrane unit



This catalogue includes information on a medical device that should be handled by qualified medical personnel and used in accordance with the instructions or label.



## Which patients can you recommend **AQUAI** to?

### RELAXATION

Modern pace of life, work and family problems make many of your patients more stressed than ever. They take too little time for themselves. With **AQUAI** you can offer them an ideal way to relax along with various health-boosting effects.

### PREVENTION

**AQUAI** offers a beneficial service for patients especially interested in their health and well-being, that keep their body in shape, eat healthy and try to boost their immune system and metabolism.

### THERAPY

Recent studies found that over 70% of our population has backache or other chronic conditions. They require your particular attention while **AQUAI** provides an ideal addition to your prescribed treatment.

### RELAXING MUSCLES

Highly repetitive activities, sitting, standing still, all kinds of office work or physical exertion usually lead to backache. They occur both in shoulder and neck areas as well as lumbar region. Tension or muscles hardening can be a reason of serious injuries such as slipped disc. **AQUAI** is a perfect solution for treating these problems.



## Designed for a **comfortable,** **soothing relaxation experience**

AQUAI has a wide range of uses in wellness and fitness. It perfectly suits rehabilitation centers. It is also a very good option for relaxation rooms in companies.

A patient in a rehabilitation center or a guest in a SPA center can experience the beneficial power of hydro massage without getting wet. The massage is done through a layer of specially designed membrane.

### Advantages of the treatment:

- Improvement in blood circulation and metabolism
- Increase in oxygen transport in muscles
- Improvement in movement range of extremities
- Tension relief
- Stress reduction and relaxation
- Well-being




Easy to get on & off



Comfortable design



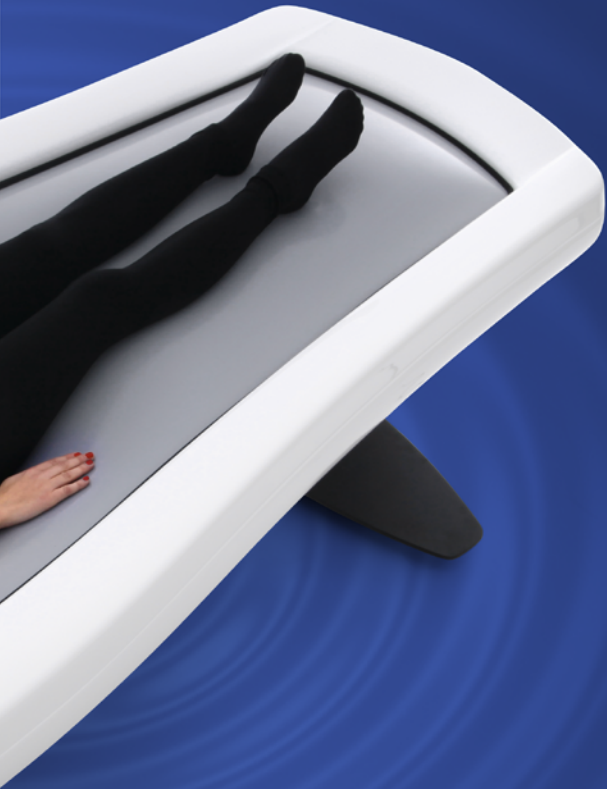

## Device features

- 6 programmable and independent zones
- All jets can work simultaneously
- Dry water massage nozzles with increased durability
- 6 standard programs
- 1 user program
- Membrane surface temperature setting adjustable from 30 to 40°C
- Adjustable treatment time (1 - 30 min)
- Delayed start at programmed preset times
- Adjustable frequency and sequence of zones
- Integrated heating system
- Internal water cooling system maintains water temperature during preset hours regardless of the number of treatments
- No need to connect to water tap permanently
- Easily accessible modern colour touch panel
- Ambient LED lighting at the base



### Comfortable water temperature

30 to 40°C



### Easy installation

230 V



### Soft music from two speakers

Audio option



### Motion sickness protection

Unique membrane

Imagine SPA interior different from the others.  
The marriage between pure beauty and ergonomomy.  
Imagine the elegance inspired by nature's harmony.  
Welcome to the **AQUAI world**.




One-touch start






## User-friendly interface

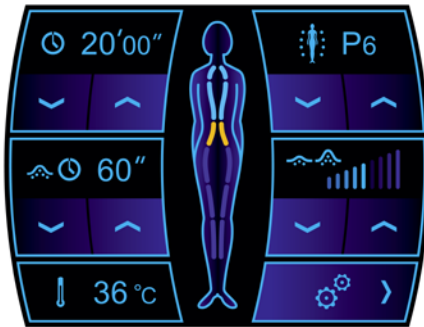


### Main screen

- Adjusting settings before and during treatment session
- Selecting required program, length of the treatment, sequence of zone changes and other settings during the session
- Welcome screen after 3 minutes of inactivity

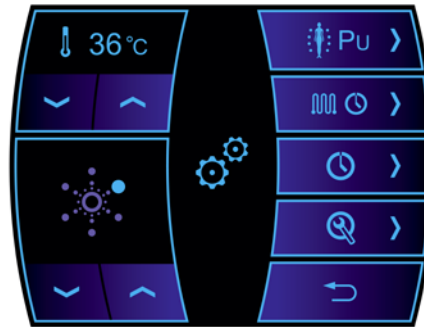

### Configuration screen


- Water temperature adjustment from 30 to 40°C
- Bottom LED light colour selection
- Easy access to User Program edit screen
- Automatic preparation setup screen
- Quick clock setup access

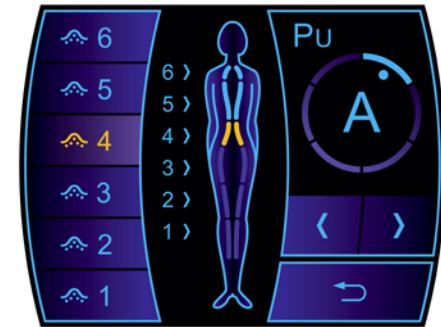

### User program edit

- Possibility of zone selection for each step
- Displaying currently working zones (yellow) and subsequent zones (blue)
- Up to 6 zones working in a sequence of a program
- Up to 3 zones working at the same time


## 6 independent zones

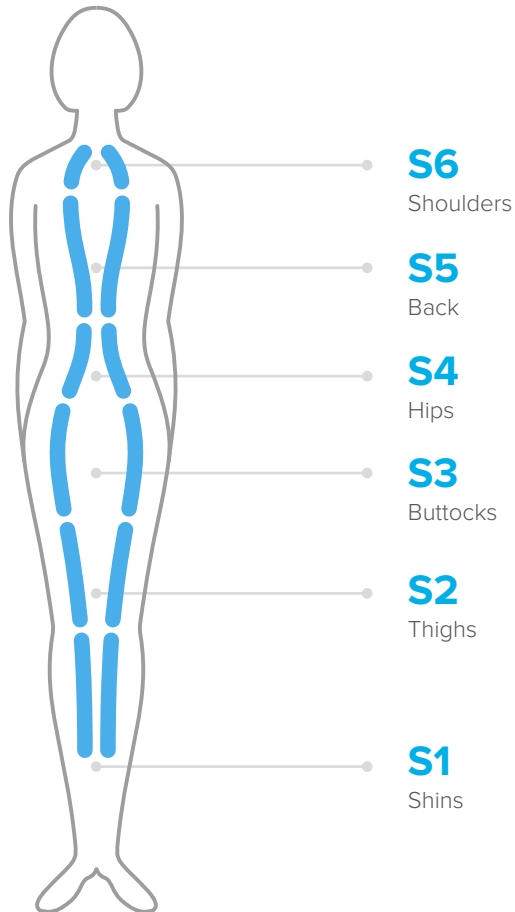

Discover the new level of dry hydromassage.


Now it is for you to decide which body part is massaged. Easily program each of the 6 fully independent working zones for better treatment results. Entire body hydromassage for deep relaxation, shoulders for neck pain relief, thighs and shins for movement improvement.

Step into the revolutionary world of AQUAI and provide your patients the greatest treatment experience ever!

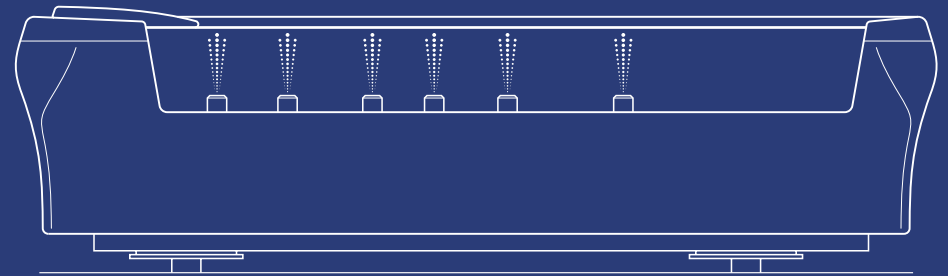



Effective & powerful  
massage

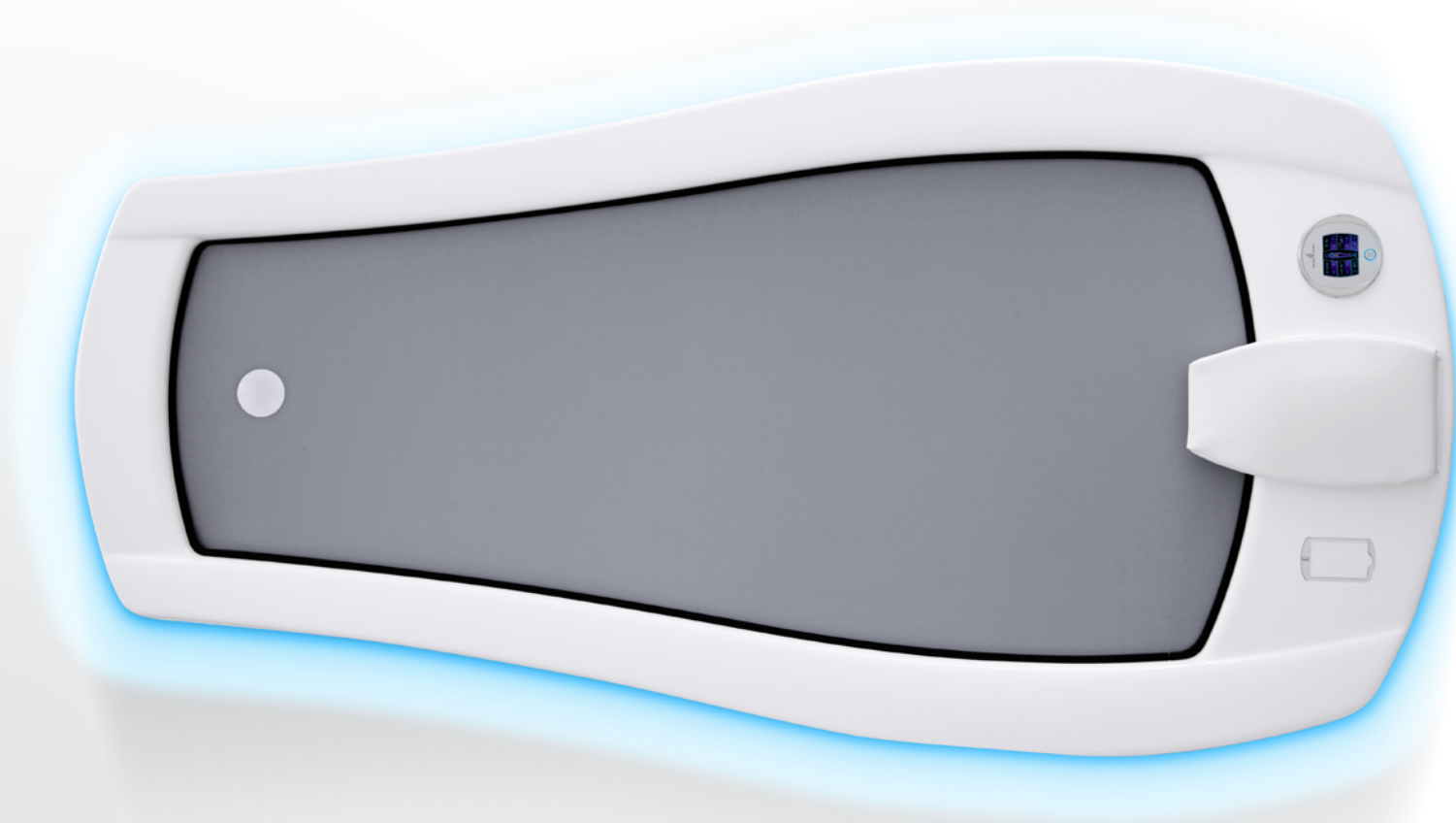

## Additional options



### CARD MANAGER


A system supporting the sale of treatment sessions using stored-value cards; together with top-up software.



### AUDIO

High-quality speaker system with sound input from an external device such as HiFi, MP3 or Bluetooth.



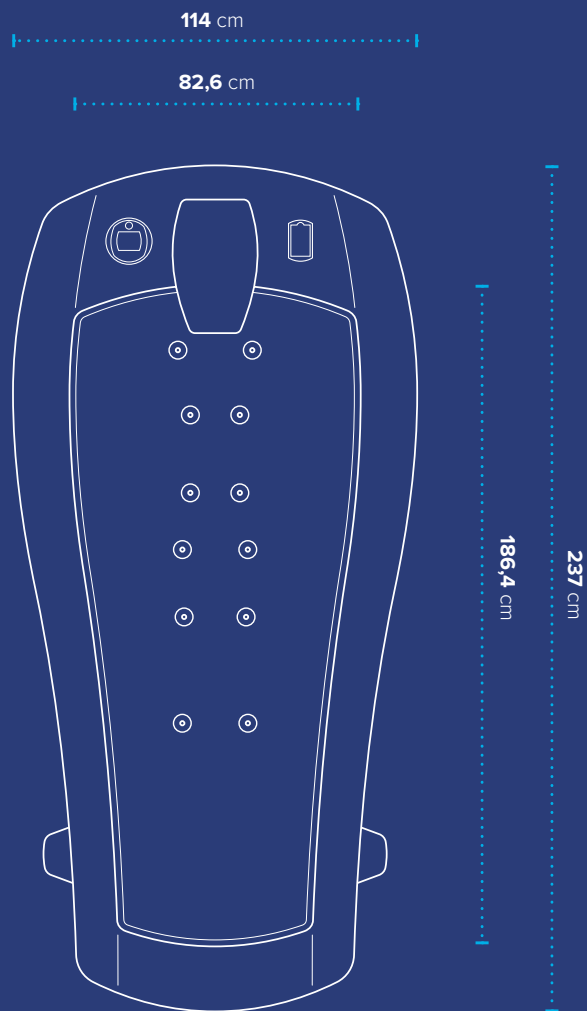

## Technical data

Dimensions (L x W x H) [cm]:	237 x 114 x 70
Capacity [l]:	220
Weight [kg]:	190
Max. load [kg]:	135
Basin colour:	White
Protection class:	I
Power supply [V/Hz]:	230/50 or 400/50

The basin and its enclosure are made of polyester resin, reinforced with glass fibre. These materials guarantee long-lasting, trouble-free operation of the device. The water system, with the exception of connection pieces and valves, is made of PVC, which ensures its reliability.



Length: **237 cm**



Height: **70 cm**



Width: **114 cm**



Max. load: **135 kg**



Weight: **190 kg**





**Meden-Inmed Sp. z o.o.**

Wenedów 2, 75-847 Koszalin, Poland

Tel.: +48 94 347 10 50 / 53

fizjoterapia@meden.com.pl

www.en.meden.com.pl

