

Neurological rehabilitation table

SAFARI ELEPHANT
B90-S1 / B90-S2 / B140-S1



This leaflet contains information on, among other things, medical devices that should be handled by qualified medical personnel and used in accordance with the instructions for use or the label.

Safari Elephant

B90-S1 / B90-S2 / B140-S1

Neurological rehabilitation table Safari type Elephant is intended for physical treatment of children with central nervous system disabilities and for performing Bobath and Vojta therapy and adults with central nervous system disorders and limited mobility. To accommodate varying needs and expectations we provide a diverse selection of models.

Structures

1-section tables



Safari Elephant B90-S1.F0



Safari Elephant B90-S1.F0
+ foot remote controller option
+ folding barriers option
+ additional accessories: paper holder



Safari Elephant B90-S1.2F2




Safari Elephant B90-S1.2F2
+ foot remote controller option
+ folding barriers option
+ additional accessories: paper holder



Safari Elephant B140-S1.F0



Safari Elephant B140-S1.F0
+ foot remote controller option
+ folding barriers option
+ additional accessories: paper holder

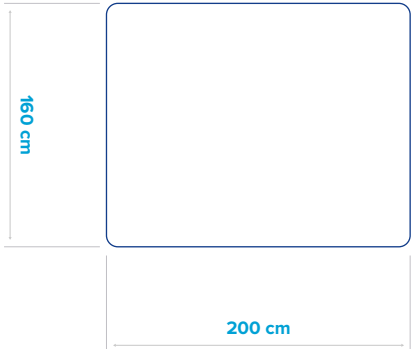




Safari Elephant B140-S1.2F2



Safari Elephant B140-S1.2F2

+ foot remote controller option
+ folding barriers option
+ additional accessories: paper holder

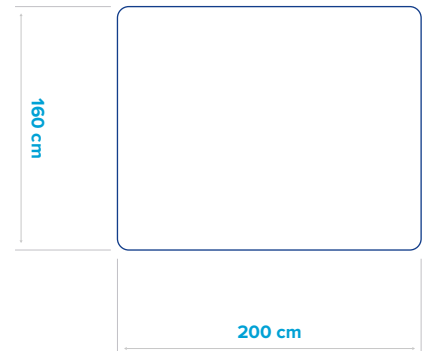

2-section tables




Safari Elephant B90-S2.F0



Safari Elephant B90-S2.F0

+ foot remote controller option
+ additional accessories: paper holder

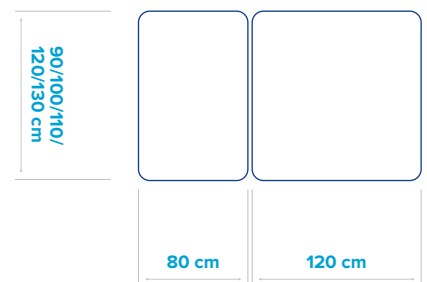


Safari Elephant B90-S2.2F2



Safari Elephant B90-S2.2F2


+ foot remote controller option
+ additional accessories: paper holder

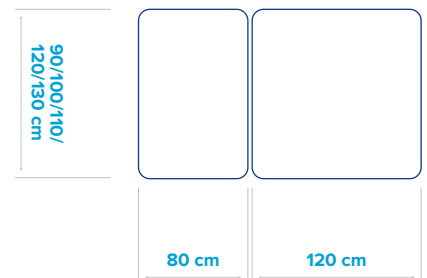

LUX upholstery:


Surface composition: 100% Polyurethane

Backing composition: 100% Polyester

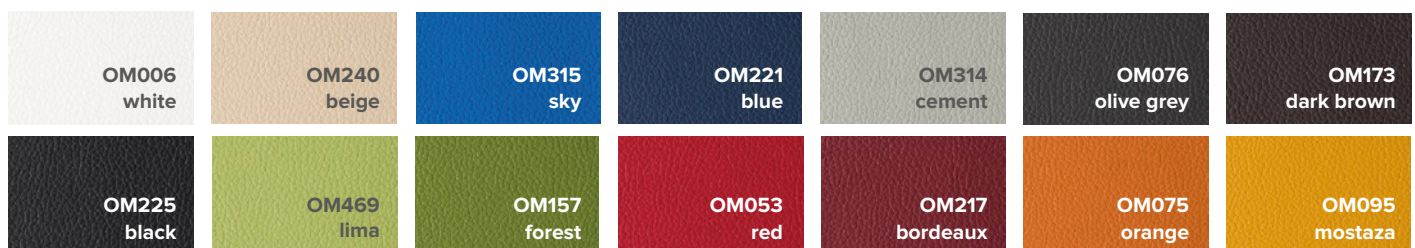

Features:

- 1- or 2-sections tables.
- Adjustable head section 0° – 85° via 2 gas springs (2-sections models B90-S2).
- A wide range of models designed to meet diverse needs and expectations.
- Extra wide surface in many possible sizes:
 - 90/100/110/120/130 cm available for tables B90-S1, B90-S2,
 - 160 cm available for tables B140-S1.
- Two-layer padding with upholstery in many colours (LUX Upholstery).
- Table top sides lined with special upholstery that is highly resistant to damage.
- Powder coated, very wide and stable steel frame.
- Electric high/low height adjustment from 50 to 99 cm via hand remote controller.
- Foot remote controller for electric high/low height adjustment (option).
- Authorization via key on hand remote controller.
- AutoLock function (option) – ensures protection against accidental changes in table position. It activates automatically and locks actuators 10 s after last use of table controls.
- F0 wheeling system – composed of 2 legs with rubber feet and 2 legs with small wheels (Ø 5 cm)
- 2F2 wheeling system – dual-lever retractable castor wheeling system composed of 2 x 2 directional castors (Ø 7,5 cm) and 4 legs with rubber feet. Dual-levers at the front and at the rear of the table operate independently.
- Rubber feet with 1 cm height adjustment allowing for the levelling of the table.
- The distance between the bottom frame and the floor: 15 cm.
- Max. lifting capacity up to 200 kg.
- Max. working load up to 275 kg.

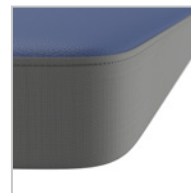

Lined with special upholstery.



Equipped as standard with a hand remote controller.



Additional accessories:

- Paper holder.
- Wedges, bolsters, semi-bolsters.
- Stools for therapist.



Technical data:

	B90-S1.F0 / B90-S1.2F2	B140-S1.F0 / B140-S1.2F2	B90-S2.F0 / B90-S2.2F2
Length [cm]:	200	200	200
Width [cm]:	90 / 100 / 110 / 120 / 130	160	90 / 100 / 110 / 120 / 130
Number of sections:	1	1	2
Height adjustment [cm]:	50-99 (electrical)	50-99 (electrical)	50-99 (electrical)
Power supply [V/Hz/A]:	100-240 / 50-60 / 3-1,5	100-240 / 50-60 / 3-1,5	100-240 / 50-60 / 3-1,5
Wheels F0 / Castors 2F2 [cm]:	5 / 7,5	5 / 7,5	5 / 7,5
Max. lifting capacity [kg]:	200	200	200
Max. working load [kg]:	275	275	275
Weight [kg]:	≤ 150 (depending on the version)	≤ 150 (depending on the version)	≤ 150 (depending on the version)



Additional options:

- AutoLock function.
- Foot remote controller.
- Folding barriers, set of 2 pcs.
(only for 1-section tables B90-S1 / B140-S1)

